


The Joerns WeCare™ bed continues our legacy of quality and reliability, now in a convenient package designed for delivery outside of institutional care. Designed and field tested for one-person delivery and setup, WeCare brings fall prevention, resident comfort and caregiver safety to any setting. Care outside the institution no longer has to be a compromise.



Enhanced Capabilities and New Options

- **UNPARALLELED STABILITY:** Enhance confidence entering / exiting the bed with full range stability and foot placement.
- **A SAFER SYSTEM:** Protect against falls injury and reduce caregiver strain with the 8.75" ultra-low height and ergonomic 26" high height. Add peace of mind knowing WeCare is guaranteed to meet bed entrapment guidelines.
- **ENHANCED COMFORT AND CARE:** Improve sleep and reduce movement in bed with a fully adjustable sleep deck and AutoContour positioning.
- **A BETTER FIT:** Improve comfort, simplify repositioning and reduce falls risk with flexible sizing. WeCare offers optional 4" length extension and width extensions from 35" to 39" (Twin XL) and 42".
- **EASIER INSTALLATION:** Protect installers from excessive strain and simplify delivery with WeCare's unique transport dolly and stair climber.
- **FLEXIBLE DELIVERY:** Receive WeCare the way you want it. Packaged for easy delivery to suit your needs, the bed system can ship palletted or by parcel carrier.

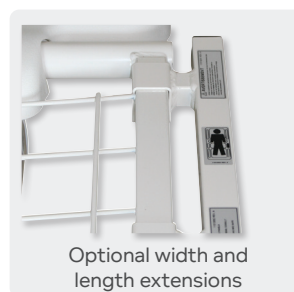

Standard Offering

- ◊ **Reduced risk of injury:** With greater travel range
- ◊ **Exceptional clearance:** Open architecture legs
- ◊ **Comfort and shear reduction:** Three-section mattress support platform
- ◊ **Better sleep deck:** Therapeutic grid surface
- ◊ **Industry leading stability:** Durable tubular steel design
- ◊ **Intuitive ergonomic hand pendant:** 8-function control with integrated hook
- ◊ **Managed controls:** Pendant function lock-out available with optional staff control
- ◊ **Safer repositioning and egress:** Entrapment compliant assist device can be mounted at head and foot
- ◊ **Effortless operation:** Four quality DC motors

Delivering care together...



We Care about delivering fast and effective care as much as you do! We've made ordering and setup of your bed quick and easy. Receive WeCare the way you want it. Packaged for easy delivery to suit your needs, with flexible delivery options. Joerns WeCare elevates patient comfort and care with simplified care solutions.

Extending a helping hand



F028 Deluxe Assist Handle



F025 Assist Handle



F026 Assist Handle



F14 Assist Device



F168 Assist Device
(requires N832AL mount)

Caring for more patients



Length Extension



Width Extension



Panel Style and Finish Options



Trapeze Patient Helper

Delivering care safely



Transport Dolly and Stair Climber

Storage and Transport Cart

Specifications

Bed Mobility	Roll-in-Low
Hand Pendant Functions	Standard: 3 Function
	Plus AutoContour
Low Height	8.75"
High Height	26"
Bed Length	Standard: 80"
	Optional: 4" Extender
Bed Width	Standard: 35"
	Optional: CareWide Accessory 39" (TwinXL) and 42"

Casters	Four casters (two locking)
Safe Working Load	450 lb
	400 lb Max Patient Weight
Underbed Clearance	18"
Product Weight	145 lb (77.7 lb foot / 67.3 head)
Electrical	120 VAC, 60 Hz, 4.0 Amps, Class II
Regulation Compliance	Conforms to IEC 60601-2-52
Warranty	1 year mechanical / electrical, 7 year structure, life welds